

Label-Free Analysis by UHPLC with Charged Aerosol Detection of Glycans Separated by Charge, Size and Isomeric Structure

*David Thomas, Ian Acworth, Bruce Bailey, Marc Plante, and Qi Zhang
Thermo Fisher Scientific, Chelmsford, MA, USA*



Overview

Purpose: To develop fast and sensitive UHPLC methods suitable to directly measure the glycan content and profile of glycoproteins.

Methods: N-linked glycans were released from proteins by PNGase-F. The native glycans were separated by ultra high performance liquid chromatography (UHPLC) on a column that employs both weak anion exchange and reversed-phase separation mechanisms, thus resolving glycans based on charge, isomerism and size. As a result, this column supports direct glycan quantification of many more isoforms than can be quantified using other column types. The native glycans were detected directly without derivatization by using charged aerosol detection (CAD).


Results: Precision, detection limits and dynamic range of quantitative measurements are presented. Figures of merit include low-nanogram (low-pmol) on-column sensitivity, over two orders of magnitude of dynamic range, and peak area precision of about two percent RSD.

Introduction

Glycoproteins of biological, diagnostic or therapeutic interest owe key aspects of their normal function to the oligosaccharides attached to the protein backbone. Changes in the number, type, composition or linkage pattern of these glycans may serve as a biomarker of disease or influence the efficacy of a biotherapeutic product.¹ For this reason, the ability to correctly identify and measure these glycans is of scientific interest, and to do so reliably, quickly and inexpensively is of practical benefit. This work explores direct detection of native glycans as an alternative to the common techniques for glycan analysis that rely on derivatization reactions to render glycans detectable. The lack of a detectable chromophore in native glycans is overcome by using HPLC with charged aerosol detection, a detector that can quantitatively measure any non-volatile compound.

Figure 1 depicts the operation of the charged aerosol detector. At the top left (1) the liquid mobile phase from the LC column enters the detector, where it is nebulized by combining with a concentric stream of nitrogen gas or air (2). The fine droplets are carried by bulk gas flow to the heated evaporation sector (3) where desolvation occurs to form dry particles (5) from any nonvolatile or semivolatile species. Any remaining large droplets drain away to waste (4). The dry analyte particles exit from evaporation and combine with another gas stream that has been charged by a high voltage Corona charger set near 2 kilovolts (6). The charged gas mixes with the dry particles and transfers positive charge to the analyte particle's surface (7). After passing through an ion trap (8) that removes any high mobility species, the remaining charged particles pass to a collector and are measured by a sensitive electrometer. The signal produced (9) is directly proportional to the quantity of analyte.

FIGURE 1. Charged aerosol detector and principle of operation.

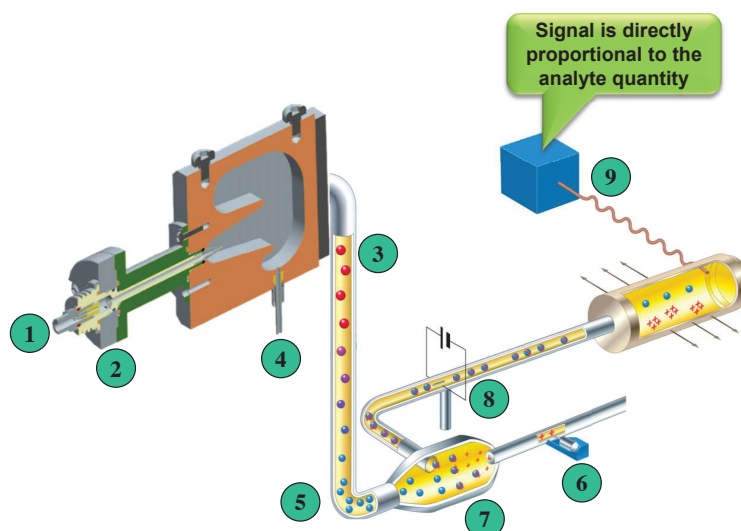

Methods

Liquid Chromatography

Thermo Scientific™ Dionex™ UltiMate™ 3000 RSLC system with:

- Thermo Scientific™ Dionex™ Corona™ Veo™ Charged Aerosol Detector:
 - Evaporation Temperature: 35 °C
 - Power function: 1.00
 - Data collection rate: 2 Hz
 - Signal Filter: 3.6 sec

Reagents: HPLC- or LCMS-grade or better

Data Analysis

Thermo Scientific™ Dionex™ Chromeleon™ Chromatography Data System (CDS) 7.2

Separation:

Column: Thermo Scientific™ GlycanPac™ AXR-1 column,
1.9 μ m, 2.1 \times 100 mm
Column Temp: 30 °C
Flow Rate: 0.4 mL/min
Injection Vol.: 2 μ L
Mobile Phase A: Deionized water
Mobile Phase B: 100 mM ammonium formate pH 4.4
Gradient: Time, %B: -5, 4; 0, 4; 25, 29 (slope = 1 mM/min).

Sample Prep:

OligoStandard;

Add 250 μ L of HPLC grade water to one vial of OligoStandard™ Sialylated Fetuin N-linked alditols (Thermo Fisher Scientific P/N 043604). Vortex to dissolve and transfer to a plastic glass HPLC autosampler vial.

Glycan Standards:

Sialylated glycan standards were purchased from Prozyme. Monosialylated A1 (GKC-124300), di-sialylated A2 (GKC-224300), and tri-sialylated A3 (GKC-335300).

Reconstitute one vial (10 μ g) of each glycan standard with 50 μ L deionized water.

Vortex to mix and transfer to plastic autosampler vial for direct injection.

Alpha acid glycoprotein (Sigma G3643) and fetuin from fetal bovine serum (ICN) were prepared by dissolving 4 mg +/- 1 mg in 1 mL HPLC grade water.

Protein PNGase F digestions were performed by using QA-Bio PNGase F Deglycosylation kit (QA-Bio P/N E-PNG01) per the manufacturer's instructions. Briefly, add 35 μ L of protein solution to a plastic centrifuge tube. Add 10 μ L 5x Reaction Buffer 7.5 and 2.5 μ L of Denaturation Solution. Heat at 100 °C for 5 minutes. Cool. Add 2.5 μ L of Triton X-100 and mix. Add 2.0 μ L of PNGase F to the reaction. Incubate 18 hours at 37 °C. Centrifuge at 6720 x g for 10 min and inject the supernatant.

Method Development

Three method parameters were optimized during development of this UHPLC-CAD method. First, from a starting concentration of 4 mM ammonium formate, the gradient slope was optimized by comparing glycan resolution and total run time for gradient slopes ranging from 0.5 – 3 mM/min. The optimum gradient slope was 1 mM/min, as seen in Figure 2, so this was chosen for the final method.

Secondly, mobile phase composition was optimized by examining the S/N of analytes after inclusion of from 0 – 20% acetonitrile or methanol in the mobile phase. Although the signal increased with increasing organic solvent, the S/N ratio varied only slightly and peaked at 1 or 5 % (Figure 3). Given the only modest increase in S/N provided by added solvent, we chose to use a purely aqueous gradient for simplicity and robustness.

Finally, the effect of evaporation tube temperature was considered by examining S/N for the glycan analytes at evaporation temperature settings of 35, 50 and 80 °C. Although S/N clearly decreased at 80 °C, the differences between 35 °C and 50 °C were less pronounced (Figure 4). 35 °C was chosen as the evaporation temperature for the final method.

FIGURE 2. Optimizing mobile phase gradient slope.

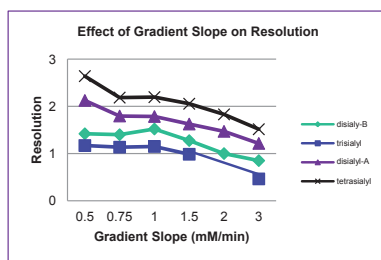


FIGURE 3. Optimizing mobile phase composition.


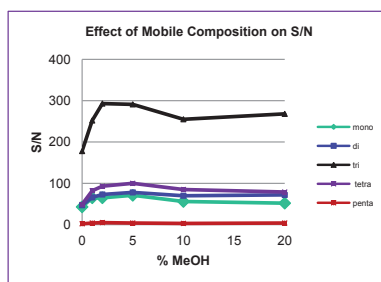


FIGURE 4. Optimizing evaporation temperature.

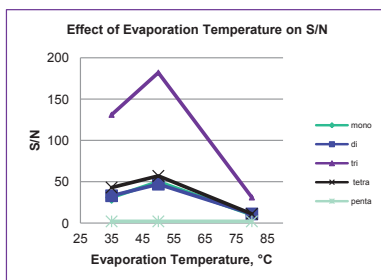


High Resolution Separation by Charge, Size and Structure

Quality control labs profile a protein's glycan pool to assess lot-to-lot variability, degradation or level of impurities. In the separation shown in Figure 5, glycans in a standard mixture are separated according to charge, size and structure. The native glycans are separated by UHPLC using a binary gradient consisting of water and a volatile ammonium formate buffer and measured directly by using a Corona Veo RS charged aerosol detector. The elution order is neutral glycans first, followed by glycans with a single negative charge (monosialylated), glycans with two negative charges (disialylated), and so on. The glycans within each charge group are separated by ion-exchange interactions. Within each well-separated charge group, glycans differing in size or isomeric structure are further resolved by reversed phase interactions.

There is no need to use fluorescent labeling when using charged aerosol detection, as may be necessary with other means of detecting these compounds. Because of the uniform response of the charged aerosol detector, the relative peak area accurately reflects the amount (pmol) within each charge group.

Note that under these conditions the neutral glycans are not well separated from the void peaks. To better resolve neutral glycans, use a shallower gradient, or derivatize the glycans to increase hydrophobicity (e.g., 2-AB)² or to introduce a negative charge (e.g., 2-AA).

FIGURE 5. Direct charged aerosol detection of bovine fetuin N-linked alditols in a standard mixture separated with high-resolution on the GlycanPac AXH-1 column.

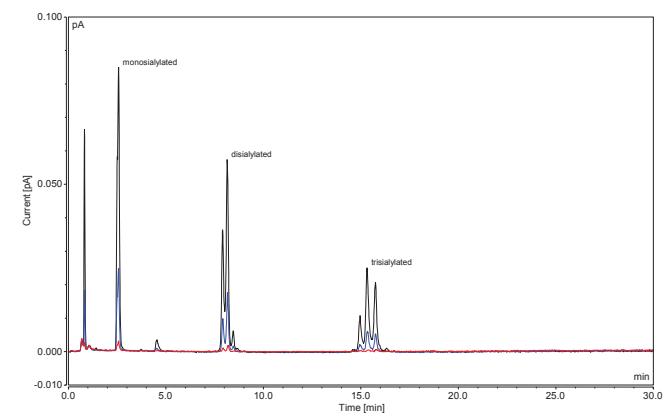

Performance

Calibration curves for three major charge groups of the bovine fetuin N-linked alditol standard (analyzed in triplicate) are presented in Figure 6. Standards were prepared at concentrations ranging from 0.5 to 100 pmol/μL; using 2 μL injections the mass on column ranged from 1 to 200 pmol. The data were fit to a quadratic equation, yielding coefficients of determination, R², greater than 0.999 for all three analytes.

Table 1 presents a summary of the method's performance, including the coefficients of determination and the limits of detection for the three major charge groups of the bovine fetuin N-linked alditol standard. Also included is precision obtained for 7 replicates of an oligosaccharide standard (110 pmol total sialylated fetuin N-linked alditols).

FIGURE 6. Calibration data for direct detection of N-linked glycans by HPLC-Charged Aerosol Detection.

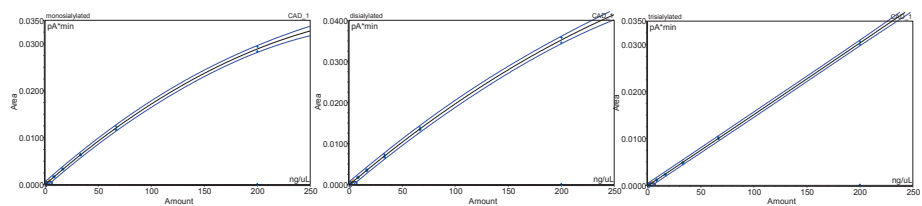


TABLE 1. Method performance for Direct Detection of Glycans by HPLC-Charged Aerosol Detection.


Component	Amount (ng/μL)	Amount (pmol)	Ret. Time ¹ (%RSD)	Peak Area ¹ (%RSD)	LOD ² (ng/μL)	R ² *
Monosialylated	4.8	5	0.22	2.35	6.5	0.9994
Disialylated	28	25	0.45	1.18	7.2	0.9993
Trisialylated	81	56	0.12	2.14	7.1	0.9996
Tetrasialylated	6.5	4	0.05	2.95	--	--
1. for n = 7 replicates						
2. Hubaux-Vos method						
* 8 levels, in triplicate (duplicate for highest concentration), quadratic fit with no offset						

Protein Digests

Charged aerosol detection is clearly highly sensitive, able to detect glycans at the low picomole concentration level. Because charged aerosol detection is universal, there might be concern that the reagents and reaction products remaining after endoglycosidase treatment would interfere with detection of the released glycans. To test for such interference, we treated two proteins by using a commercial glycan release kit and then analyzed for glycans by using UHPLC-CAD. Figure 7 shows that although additional peaks related to the reaction procedure are evident, they are well resolved from the glycan analytes and do not interfere with reliable quantification of the glycans.

This demonstrates that where ultimate sensitivity is not required³, UHPLC-CAD obviates the need to spend time or money on 2-AB labeling and delivers a clean chromatogram with no concern for reaction side products.

FIGURE 7. Direct charged aerosol detection of α_1 -acid glycoprotein and bovine fetuin N-linked glycans released by PNGase F separated with high-resolution on the GlycanPac AXH-1 column

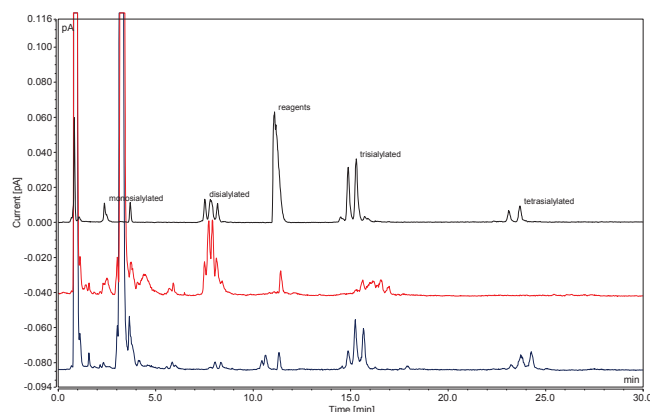

Conclusion

- The HPLC method developed to measure native glycans is precise, with retention time precision better than 0.2% RSD and peak area precision averaging 2.2 % RSD for the major sialylated N-glycans of bovine fetuin.
- Charged aerosol detection enables sensitive, direct measurement of glycans with no need to perform labeling reactions. Detection limits for native glycans are in the low pmol (ng on-column) range.
- By responding directly to any non-volatile compound, charged aerosol detection is able to measure both native and fluorescently labeled glycans, yielding simple, accurate and precise estimates of relative concentration even in the absence of pure primary standards.

References

1. Varki, A. Biological Roles of Oligosaccharides: All of the Theories Are Correct. *Glycobiology*, **1993**, *3*, 97-130.
2. Aich, U.; Saba, J.; Liu, X.; Thayer, J.; Pohl, C. Separation of 2AB-labeled N-linked glycans from bovine fetuin on a novel ultra high resolution mixed-mode column, **2014**, Thermo Scientific Applications Note 20908.
3. Hurum, D.; Basumallick, L.; Rohrer, J. Evaluating protein glycosylation in limited-quantity samples by HPAE-PAD, **2013**, Thermo Scientific Application Note 1050.

www.thermoscientific.com/chromatography

©2015 Thermo Fisher Scientific Inc. All rights reserved. ISO is a trademark of the International Standards Organization. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries. This information is presented as an example of the capabilities of Thermo Fisher Scientific products. It is not intended to encourage use of these products in any manners that might infringe the intellectual property rights of others. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.



Thermo Fisher Scientific,
Sunnyvale, CA USA is
ISO 9001 Certified.

Africa +43 1 333 50 34 0
Australia +61 3 9757 4300
Austria +43 810 282 206
Belgium +32 53 73 42 41
Brazil +55 11 3731 5140
Canada +1 800 530 8447
China 800 810 5118 (free call domestic)
400 650 5118

Denmark +45 70 23 62 60
Europe-Other +43 1 333 50 34 0
Finland +358 9 3291 0200
France +33 1 60 92 48 00
Germany +49 6103 408 1014
India +91 22 6742 9494
Italy +39 02 950 591

Japan +81 6 6885 1213
Korea +82 2 3420 8600
Latin America +1 561 688 8700
Middle East +43 1 333 50 34 0
Netherlands +31 76 579 55 55
New Zealand +64 9 980 6700
Norway +46 8 556 468 00

Russia/CIS +43 1 333 50 34 0
Singapore +65 6289 1190
Sweden +46 8 556 468 00
Switzerland +41 61 716 77 00
Taiwan +886 2 8751 6655
UK/Ireland +44 1442 233555
USA +1 800 532 4752

Thermo
SCIENTIFIC

A Thermo Fisher Scientific Brand